



**SWING**

*Simple bathing while  
sitting and lying*






## SWING


### Safety and comfort for people with limited mobility

The SWING sitting/lying bathtub is an excellent alternative to conventional bathing systems; it enables entry from the side and possesses electric lowering and tilting apparatus. The bathtub offers safety and comfort when bathing to those needing care and supports staff during daily care.

The SWING bathtub is characterized by gentle, formfitting contours with comfortable head support. In the tilted position, the patient is optimally supported and the caretaker can look after the patient without unnecessary strain. At the push of a button, the bathtub tilts to the desired working height, thus preventing additional back strain. The air jet system is stimulating and relaxing. The gears for the lifter are situated in the concavity of the base, thus allowing the patient to be optimally positioned in the bathtub.

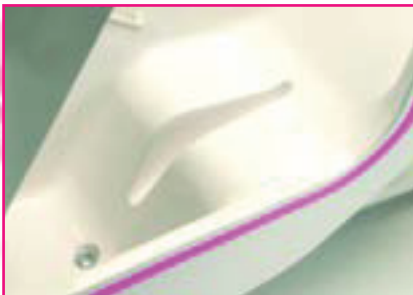

The sitting/lying bathtub is fitted with a thermostatic mixer with high temperature safeguard and corresponds to the most recent international guidelines. Additionally, a shower head with adjustable water stream is fixed to the rim of the bathtub. Those needing care can be placed directly from a wheelchair into the bathtub with a lifter. For a full bath, the water is poured in with the door open and flows from the foot basin, which fills up first, around the end of the bathtub. After closing the door, the bathtub can be tilted to a lying position. Alternatively, the bathtub can be used as a comfortable, easily accessible shower and the foot basin likewise used for foot care.

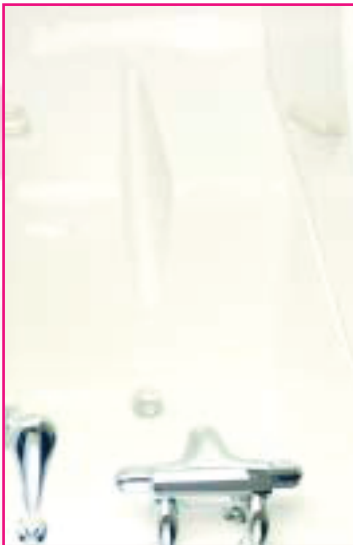

In order to avoid cross infections, it is recommended that the bathtub is equipped with an integrated disinfection system (for bathtub and whirlpool bath).

### Additional fittings

Massage bath with 10 jets and heating

Disinfecting system

Automatic bathtub filling

Your distributor:



Barrier Free Lifts™  
1620 SW 17th St  
Ocala, Florida 34471  
USA

Phone: 800-582-8732  
Fax: 866-378-3318  
[www.barrierfreelifts.com](http://www.barrierfreelifts.com)  
[info@barrierfreelifts.com](mailto:info@barrierfreelifts.com)