

# Nursing bathtub iTub SWING



Simple bathing while sitting and lying

Safety and comfort for people with limited mobility

The SWING sitting/lying bathtub is an excellent alternative to conventional bathing systems; it enables entry from the side and possesses electric lowering and tilting apparatus. The bathtub offers safety and comfort when bathing to those needing care and supports staff during daily care.


The SWING bathtub is characterized by gentle, formfitting contours with comfortable head support. In the tilted position, the patient is optimally supported and the caretaker can look after the patient without unnecessary strain. At the push of a button, the bathtub tilts to the desired working height, thus preventing additional back strain. The air jet system is stimulating and relaxing. The gears for the lifter are situated in the concavity of the base, thus allowing the patient to be optimally positioned in the bathtub.

## Nursing bathtub iTub™ SWING

The innovative safety thermostat control for bathtubs

- 3 Temperature Indicators for
  - temperature in the bathtub
  - temperature in the bathtub inflow
  - temperature in the hand shower

All functions integrated in one control panel

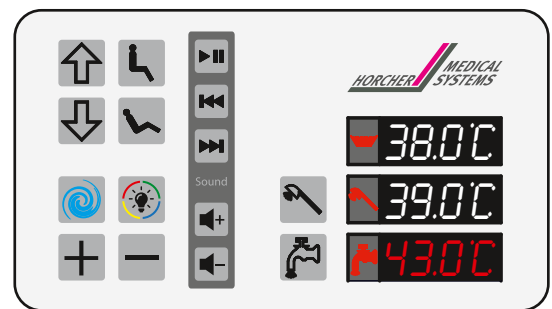

Optional

- AirSpa - 8 hydromassage airmix jets
- SoundSystem - Bluetooth or MP3 via USB
- LightSystem - 2 RGB LED spots
- AutoFill - Automatic bathtub filling
- Disinfection system
- Automatic height adjustment

Length sitting position	2040 mm
Length lying position	1960 mm
Product width	800/740 mm
Water pressure	2 - 6 bar
Water capacity	approx. 100 liter
Max. patient weight	150 kg / 330 lbs.

Barrier Free Lifts, Inc. | 1620 SW 17th St | Ocala, FL 34471

Phone: 1-800-582-8732 | Fax: 1-866-378-3318 | us-office@barrierfreelifts.com | www.barrierfreelifts.com

**BARRIER  
FREE LIFTS**